# SafetyWasteCap 2.0 S 60/61 Thread

# Universal




**THE System** for Solvent Safety

### **ORIGINAL SAFETY:**

### **Innovative & universally usable**

- SafetyWasteCap with S 60/61 thread in the new design of generation 2.0.
- · More connections thanks to the beveled core - for more applications!
- · Maximum flexibility and equipment on a single cap - no adapter required!
- · Connections which are not required can be closed with blind plugs.
- Large scope of delivery: 9x Fittings: 3x 1.6 mm, 3x 2.3 mm, 3x 3.2 mm 3x Blind plugs for capillary connection 3x Tube connectors 3x Blind plugs for tube connections

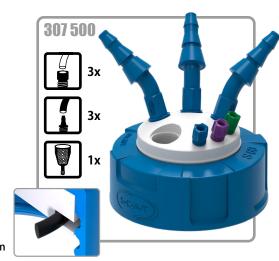

# **NEW: Ergonomic Design**

Enhanced stability and construction / Improved handling / Easy installation / Smart exchange

### **NEW: Medical Grade PTFE**

• Excellent chemical resistance / Food safe / FDA compliant / Flammability classification V-0 (UL-94)

# **NEW: PPS Screw Cap**

• Autoclavable / sterilisable up to 200 °C / Flammability classification V-0 (UL-94)

# **Proven: PTFE Body with integrated FKM Seal**

• Provides hermetical connection to the container opening and adapts smoothly even to rugged glass surfaces. An additional PTFE layer protects the seal from direct contact to solvents.





Phone: +49 (0) 61 05 / 30 55 86 - 0 **₹** e-Mail: info@scat-europe.com **♦ Internet:** www.scat-europe.com

> SCAT Europe GmbH Waldecker Straße 7 64546 Mörfelden-Walldorf Germany







Part No.	Description	Thread size	Connections for Ø 1.6 / 2.3 / 3.2 mm OD tubing	Connections for Ø 6.4 - 9 mm ID tubing	Service life (in months)	Unit
307 500	SafetyWasteCap S 60/61, V2.0, Universal	S 60/61	3x	3x	-	1
490 336	Exhaust filter size M, V3.0, with splash protection and change indicator <b>Economy package</b>	GL 14	-	-	2x 6	2